



HYPNEA 3

CIDELEC polygraph

Listening to
understand better


CIDELEC

Explore the unique functions of Hypnea 3 and how they are redefining sleep studies.

Unique

Associated with PneaVoX®, Hypnea 3 integrates the raw re-listenable tracheal sounds for improved understanding of the pathology and easier patient education.

Intuitive

A solution designed for the patient. Its intuitive interface and tutorial video adapt to all profiles for fully autonomous placement.

Modular

Thanks to an integrated Bluetooth interface, Hypnea 3 can acquire signals from optional modules (Leg EMG, machine recording, etc.)

Technological

Perfect combination of design and cutting edge technology, Hypnea 3 integrates the HypnoLight function.

Compact


Its compact format, lightness and silicone wristband ensure optimum comfort, guaranteeing reliable data all night long.

Durable

Easy to repair, requiring less energy and equipped with reusable sensors, Hypnea 3 minimizes its environmental impact.



Technical specifications

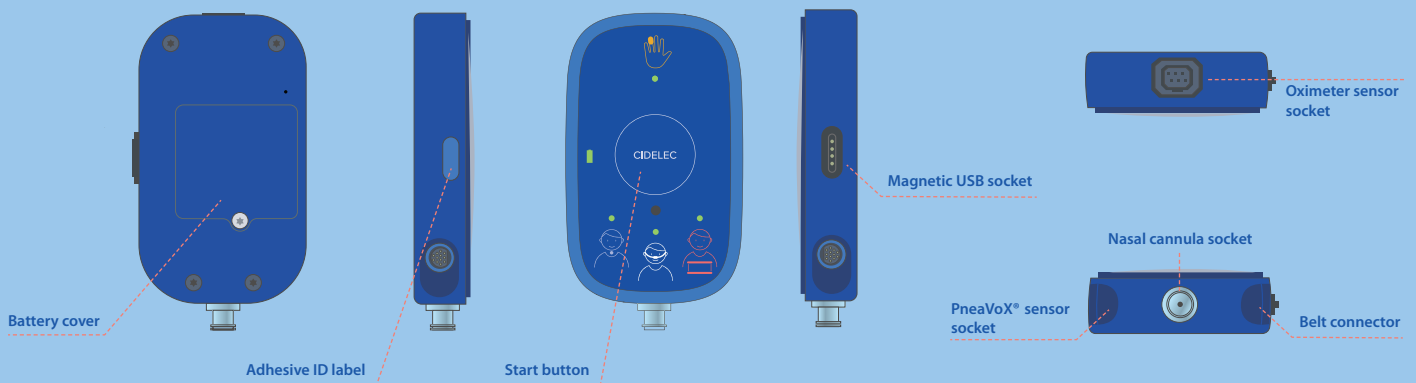

Specifications

Signals availables	Suprasternal pressure, breathing flow, nasal pressure, thoracic and abdominal efforts, snoring (sic), raw sounds, brightness, position, actimetry, SpO2, BPM, PPG.
	Optional: HLT ExG, EMG, PTA, Thermocouple
Raw sound signal	2KHz
Wireless interface	Bluetooth® V5.0 (BLE)
Memory capacity	256 Mo, transfert time 20 secs (excluding raw sounds)
Power supply	720 mAh LiPo battery, charges in 4 hours

Our works and scientific publications:



PneaVoX® :
www.cidelec.net/pneavox/



HypnoLight :
www.cidelec.net/hypnolight/



Charge Hypoxique :
www.cidelec.net/charge-hypoxique/

www.cidelec.net